

Eye Wash Station, WINLAB®



WW5PT021

WINLAB®


Eye wash station

Very easy to use eye wash station in the event of an emergency.

When the eye wash station is opened via the large arm lever, flushing begins automatically, thus freeing the individual's hands so they can hold their eyelids open, for example. The eye wash station must be connected to a cold water line that supplies drinking water. The power of the water stream can be set via a ball valve when mounting.

Features:

- Extremely wide spray heads that not only cover both eyes, but the entire face (even those models with only one spray head) with a gentle water stream
- Supple rubber prevents eyes from being further injured if they come into contact with the molding
- Spray heads always remain clean, as cap seals automatically pop off during the wash cycle
- Station meets DIN EN 15154-2:2006

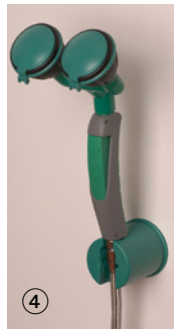

Article	Order No.
① Compact eye wash station 1/2 IG	10 03 00084
For connecting to existing outlets. With integrated BALLOFIX ballcock for regulating the jet, plus self-adhesive notice. With 2 shower heads (2 nozzles).	

Eye wash bottles

Manual eye wash stations with locking release lever and shower head. Due to one-hand operation, one hand is free in case of an emergency.

With 1.5 m silicon hose covered with stainless steel, integrated flow rate limitation and backflow prevention. Thanks to the 1.5 m tube, the hose can also be used on a reclined person.

Article	Order No.
② Manual eye wash, with 1 spray head, sprays vertically upwards	36 43 00080
Version for wall fixing and table fixing	
③ Manual eye wash, with 1 spray head, angled by 45°	36 43 00062
Version for wall fixing and table fixing	
④ Manual eye wash, with 2 spray heads, angled by 45°	10 03 00065
Version for wall fixing and table fixing	

